

Kinéris®



Braden scale score



INDICATIONS

HELP IN PREVENTING PRESSURE ULCERS FOR PATIENT AT VERY HIGH RISK OF PRESSURE ULCER DEVELOPMENT CONFINED TO WHEELCHAIR BETWEEN 8 AND 10 HOURS A DAY

WITH OR WITHOUT SEVERE STABILITY DISORDERS DEPENDING ON MODELS ACTIVE OR NON ACTIVE PATIENT DEPENDING ON HEIGHTS

PATIENT WEIGHT: 40 TO 150 KG/ 6.2 TO 23.6 ST

CNEDIMTS (1) RECOMMENDATIONS IN ITS OPINION DATED 22/12/2009

(1) NATIONAL COMMISSION FOR ASSESSING MEDICAL DEVICES AND HEALTH TECHNOLOGIES

CATEGORY OF INDICATIONS	CUSHION TYPE
Category 3a: preventive cushion for patients with a high to very high risk of developing pressure ulcers (subject to clinical judgement and assessment scales) and/or specific risk of ischial pressure ulcers, WITHOUT asymmetrical posture	With non-motorised air and individually deformable telescopic cells, PROFIL'AIR® SINGLE compartment version
Category 3b: preventive cushion for patients with a high to very high risk of developing pressure ulcers (subject to clinical judgement and assessment scales) and specific risk of ischial pressure ulcers, WITH asymmetrical posture	With motorised air and individually deformable telescopic cells, PROFIL'AIR® DUAL compartment version

- ▶ Reinforced prevention of ischial pressure ulcers risks
- ▶ The mobility of the cells allows reducing friction and shearing forces
- ▶ Dual or single compartment for users with or without stability disorders
- ▶ The height of 7 cm is relevant for active disabled person (transfer aid)
- ▶ Removable, stretchable and breathable cover for a thermoregulation when prolonged sitting
- ▶ Inner cover with non-slip base to ensure the cushion keeps in place
- ▶ Delivered with a manometer to facilitate the inflation
- ▶ Various sizes available to meet different needs

KEY FEATURES


- ▶ Cushions made up of adjustable and individually deformable pneumatic air cells
- ▶ Dual or single compartment for users with or without stability disorders
- ▶ Several sizes and two heights available: 7 or 10 cm high
- ▶ Delivered with a manometer to facilitate the inflation
- ▶ Integral, removable cover, SMB cord in PES FR (Polyester / Polyurethane + Elastane on the surface) :
 - Breathable.
 - Stretchable.
 - Highly resistant to abrasion.
 - Resistant to mold and light.
 - Crease-resistant.
 - Anti-allergenic (Oeko-Tex® class 1).
 - Flame Retardant.
 - Washable at 40°C (Occasional washes possible at 90°C).



SINGLE-COMPARTMENT



Single-compartment for users with no stability disorders.

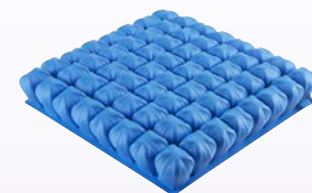


Single-compartment 7 cm.

DUAL-COMPARTMENT



Dual-compartment for users with stability disorders

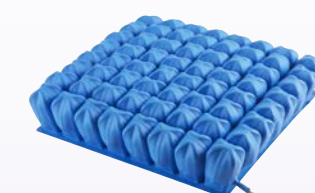

Dual-compartment 10 cm.

* The prescriber must ensure that patients or their caregivers are able to adjust the support. We recommend a visit of the prescriber the following month.

KINÉRIS® Integral removable cover Non-slip base	CUSHION HEIGHT 7 CM		CUSHION HEIGHT 10 CM		PATIENT WEIGHT KG	DIMENSIONS A X B
	SINGLE-COMPARTMENT	DUAL-COMPARTMENT	SINGLE-COMPARTMENT	DUAL-COMPARTMENT		
V367-1EL	V367-2EL	V3610-1EL	V3610-2EL	40 to 55 kg (6.2 to 8.6 st)	28.5 x 28.5 cm / 11 x 11 inches	
V427-1EL	V427-2EL	V4210-1EL	V4210-2EL	40 to 45 kg (6.2 to 7.08 st)	28.5 x 33 cm / 11 x 12.9 inches	
V497-1EL	V497-2EL	V4910-1EL	V4910-2EL	50 to 60 kg (7.8 to 9.4 st)	33 x 33 cm / 12.9 x 12.9 inches	
V567-1EL	V567-2EL	V5610-1EL	V5610-2EL	40 to 65 kg (6.2 to 10.2 st)	33 x 37.5 cm / 12.9 x 14.5 inches	
V667-1EL	-	V6610-1EL	-		36 x 36 cm / 14.1 x 14.1 inches	
V647-1EL	V647-2EL	V6410-1EL	V6410-2EL	50 to 75 kg (7.8 to 11.8 st)	37.5 x 37.5 cm / 14.5 x 14.5 inches	
V727-1EL	V727-2EL	V7210-1EL	V7210-2EL	50 to 85 kg (7.8 to 11.8 st)	37.5 x 42 cm / 14.5 x 16.5 inches	
V817-1EL	V817-2EL	V8110-1EL	V8110-2EL	55 to 90 kg (8.6 to 14.1 st)	42 x 42 cm / 16.5 x 16.5 inches	
V907-1EL	V907-2EL	V9010-1EL	V9010-2EL	85 to 120* kg (13.3 to 18.8 st)	42 x 46.5 cm / 16.5 x 18.1 inches	
V677-1EL	-	V6710-1EL	-		43 x 36 cm / 16.9 x 14.1 inches	
V777-1EL	-	V7710-1EL	-		43 x 42 cm / 16.9 x 16.5 inches	
V107-1EL	V107-2EL	V1010-1EL	V1010-2EL	110 to 120* kg (17.3 to 18.8 st)	46.5 x 46.5 cm / 18.1 x 18.1 inches	

**Compatibility subject to medical opinion up to a weight of 150 kg.

3 years warranty

Cleaning instructions:



Photos not contractually binding

FP000045/01_2018

